


Ultra Mirage™ Full Face Mask



The ultimate solution for mouth leak

Mouth breathing and leak are common problems for CPAP and bilevel therapy users. These problems result in ineffective therapy, discomfort and disrupted sleep—all of which can cause patients to be noncompliant or to abandon therapy altogether.

ResMed's Ultra Mirage™ Full Face Mask covers both the nose and mouth so patients receive effective therapy even in the presence of mouth breathing and leak.

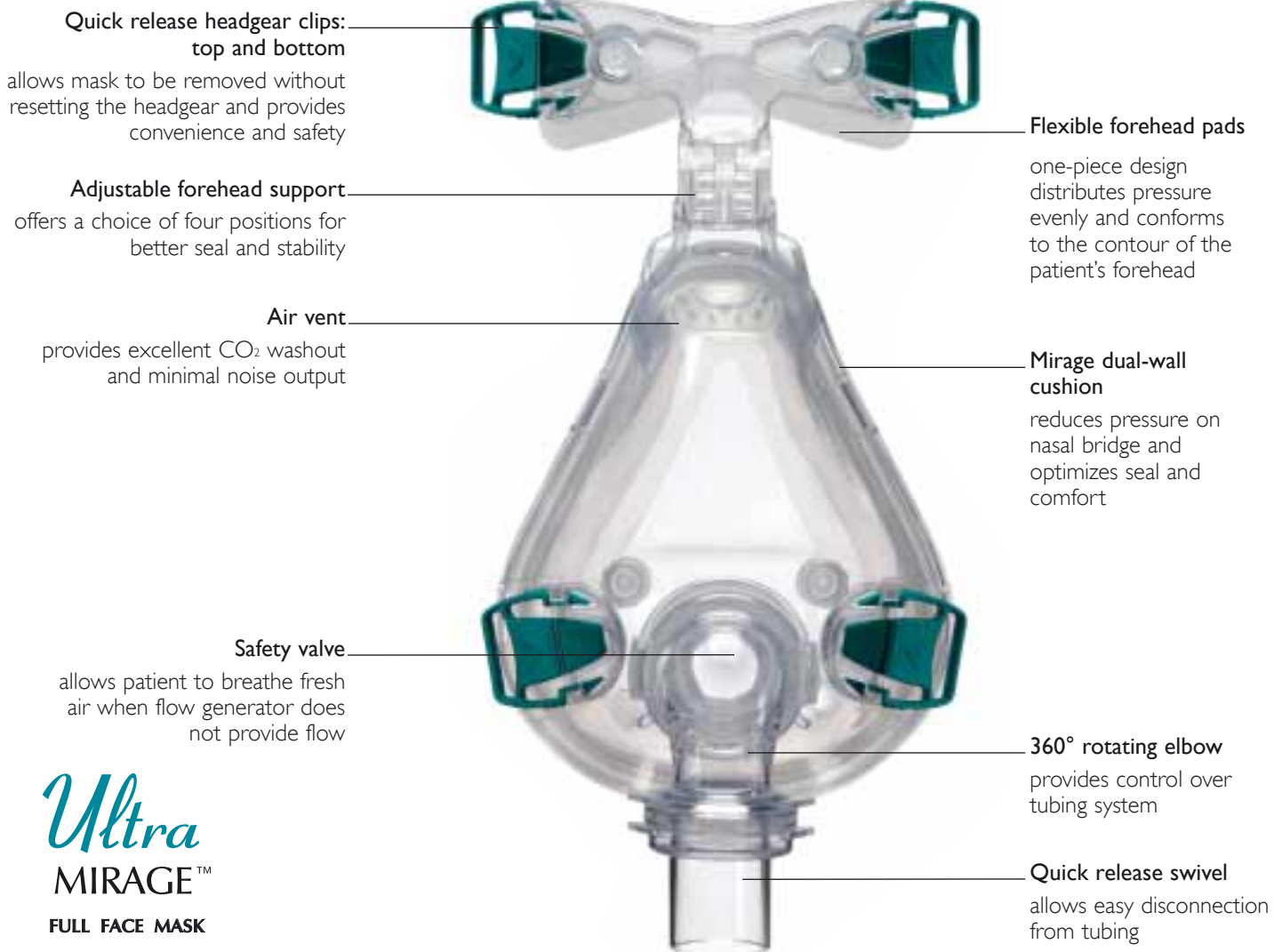

Ultra
MIRAGE™
FULL FACE MASK

Technology highlight—adjustable forehead support and pads

The adjustable forehead support takes seal performance to a new level.

- Enhanced stability is achieved by minimizing movement of the mask relative to the face. Seal is maintained even when the patient moves during sleep.
- Four adjustable positions optimize seal and comfort.
- The flexible forehead pads conform to the contour of the patient's forehead.
- The pads distribute pressure evenly over a larger surface area, reducing pressure on the patient's face.



Patient benefits

- **Trusted comfort with proven technology**
Mirage cushion technology reduces pressure on the nose and optimizes seal and comfort.
- **Superior seal through enhanced stability**
Adjustable forehead support and pads help maintain seal even with patient movement.
- **User convenience**
Quick release headgear clips enable quick and easy mask fitting and removal. Once set-up, the headgear doesn't have to be readjusted.
- **Peace of mind**
 - Minimal dead space, superior airflow characteristics, and the built-in air vent minimize CO₂ rebreathing.
 - Quick release headgear clips allow easy removal.
 - Built-in safety valve allows patient to breathe fresh air when the flow generator provides no airflow.
- **Efficient mask size selection**
 - 3-D fitting templates help achieve first-time fit.
 - Templates may also help reduce mask fitting time and minimize call backs and returns.

Professional benefits

- **Confidence in effective therapy delivery**
 - As the mask covers both the nose and mouth, it provides effective therapy in the presence of mouth breathing, and
 - Minimizes the effects of mouth breathing to help optimize patient compliance.

Ultra Mirage Full Face Mask System

Product codes

North & Latin America

Small Standard	60600
Small Shallow	60601
Medium Standard	60602
Medium Shallow	60603
Large Standard	60604
Large Shallow	60605

Europe and Asia Pacific

Small Standard	60610
Small Shallow	60611
Medium Standard	60612
Medium Shallow	60613
Large Standard	60614
Large Shallow	60615